# Bitvae

## UK REP

APEX CE SPECIALISTS LIMITED  
89 Princess Street, Manchester, M1 4HT, UK

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
Thank you for choosing the D2 sonic toothbrush. Please read the instructions carefully before using the product, and keep the user manual in a safe place for future reference.

## 1. Features:

- A Waterproof and non-slip design.
- Replaceable brush heads with unique 3D shape bristles design and a bristle replacement reminder, to help keep your teeth whiter and healthier.
- The high-frequency vibration effectively removes plaque, giving your teeth and gums a deep cleaning, and preventing gingivitis.
- USB charging, travel friendly and compatible with global voltage.
- A low-power reminding function, when the voltage of the battery is low on power, the indicator light will keep flashing quickly or the toothbrush is unable to turn on during the operating process, as a reminder to charge the toothbrush.
- A 2 minute smart timer with interval pausing, that gives a reminder to change the brushing location every 30 seconds.
- There are 5 modes: White, Clean, Soft, Polish and Gum Care settings built into the toothbrush. There is also a mode memory function built into the brush, that means when you turn on the toothbrush, the default starting mode stays the same as the mode that was set at the time the toothbrush was last shutoff.

## 2. Designs:

- ① Brush Head
- ② ON/OFF Button
- ③ White
- ④ Clean
- ⑤ Soft
- ⑥ Polish
- ⑦ Gum Care
- ⑧ Charging Indicator
- ⑨ Charging Port
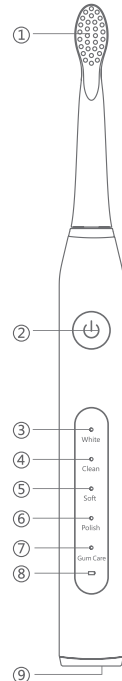
### 3. Operations:

#### ★ Using the toothbrush

- Slide the toothbrush head down to the handle and make sure it is locked in the correct position.
- Turn on the toothbrush by pressing the “Power” button.
- Shift between the 5 modes (“White, Clean, Soft, Polish and Gum Care”) by continuously pressing the “Power” button.
- Wet the bristles and apply an appropriate amount of toothpaste onto the brush.
- Place the bristles with toothpaste against the teeth at a slight angle towards the gum line.
- Gently move the brush head across the teeth on the front, back, and all chewing surfaces. The brush should be allowed to glide across the teeth with minimal pressure. Do not apply too much pressure to the teeth/gums, or apply the bristles to one area of the mouth for an extended period of time.
- The brush has a 2 minute smart timer as this is recommended by dentists as the optimal time for brushing the teeth and there is interval pausing after every 30 seconds as a reminder to move to another quadrant of the mouth.
- Move to the next section when you feel the pause and the toothbrush will automatically power off after 2 minutes.
- To avoid splashing toothpaste, only turn on the toothbrush after you have placed it against your teeth and only remove it from the mouth after it is turned off.

#### ★ Changing the toothbrush head

- Install: Ensure the point on the brush head rod beside the arrow and the toothbrush buttons facing you in the same direction. Align the point with the groove on the brush handle. Press down until it clicks to the endpoint, then rotate clockwise to fasten it.
- Detach: Hold the brush handle firmly, rotate the toothbrush head counterclockwise until it is at the endpoint, then pull it out from the brush handle.

#### ★ Charging the toothbrush

- If the toothbrush is in low battery condition when turning it on for the first time, please ensure it is fully charged before usage.
- If the toothbrush is not in use for more than 3 months, please fully charge it according to the charging instructions before using it.
- If the indicator light flashes quickly during the operation process, it indicates the toothbrush is in a low battery condition, please charge the toothbrush as soon as possible.
- When you are going to charge your toothbrush, please insert the USB cable into the charging port in the bottom of the toothbrush and ensure the brush is connected to the electric power safely.
- The charging indicator light will be white and flashes slowly during the charging process. When the battery is fully charged, the indicator will turn on for two minutes and not flash anymore, after 2 minutes, it will turn off automatically. After the toothbrush is fully charged, it is possible to use it for two minutes per session, twice a day for 60 days.
- There is a waterproof function built into the toothbrush, so it can be used in the bathroom safely.

#### ★ Modes Introduction

##### • The first mode: White

In this mode, the white indicator light turns on in the “White” mode position, the swing of the brush is large and the vibration intensity is the strongest. It is suitable for users

who have high requirements for cleaning effect and do not have sensitive teeth.

##### • The second mode: Clean

In this mode, the white indicator light turns on in the “Clean” mode position, the swing of the brush is normal and the vibration intensity is in medium speed. It is suitable for normal oral conditions, to gain the strongest cleaning effect, and maintain healthy oral hygiene.

##### • The third mode: Soft

In this mode, the white indicator light turns on in the “Soft” mode position, the brush swing is soft and the vibration intensity is the softest of the five modes, which is suitable for users with sensitive teeth.

##### • The fourth mode: Polish

In this mode, the white indicator light turns on in the “Polish” mode position, the bristles’ swing and the vibration intensity changes rapidly. The simulated effect can be used for additional treatment of the incisor area which has a whitening effect.

##### • The fifth mode: Gum Care

In this mode, the white indicator light turns on in the “Gum Care” mode position, the bristle swing and the vibration intensity changes slowly. This is for gentle stimulation, designed to keep gums firm and healthy.

### 4. Specifications:

Working temperature: 0-50°C

Working voltage: DC3.7V

Bristles material: Nylon

Input voltage: DC5.0V

Rated power: 1.5W

Waterproof rating: IPX7

### 5. Precautions:

- Do not disassemble, repair, or alter the product without authorization.
- People who have severe periodontal disease or oral surgery should consult a doctor before use.
- Toddlers, infants and people who are unable to operate should be supervised by a parent or adult.
- Do not soak the product in water for an extended period of time.
- Do not place the toothbrush in high temperature condition or expose it to the sun.
- Do not force the bristles hard against the teeth or gums.
- Not suitable for children under 12 years old.

### 6. After-sales service:

- A 1-year warranty starts on the day of purchase, please contact us if it is necessary.
- Within the warranty period, any defect arising from materials or a product malfunctioning, we will replace the parts or the whole product depending on the situation and its condition.
- When you apply for a warranty based repair, you must submit the toothbrush and proof of the date of purchase (purchasing invoice) to us.
- The following conditions are not included in the warranty:
  - Damages caused by improper use, for example, improper use of the power supply, plugged into improper socket, artificial damage, etc.
  - Ageing and wearing-out caused by normal use.
  - Damages due to product being used for other purposes, or use of components not supplied by us.

There is a built-in rechargeable lithium battery in the product. When disposing of the appliance, do not mix it with general household waste. It should be sent to the battery recycling center for disposal in accordance with local disposal management regulations.

### **Recycling & Environmentally Friendly Disposal**

In order to protect the environment, after the expiration of the products service life, it should not be discarded in the same way as domestic garbage; please recycle at a professional recycling center. For further information on the recycling of this product, please consult your local municipality, processing service center, or the store where you purchased the product. You can help protect the environment! Please remember to respect the local regulations, and hand in the non-working electrical equipment to an appropriate waste disposal center.

## **IMPORTANT SAFEGUARDS**

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

**READ ALL INSTRUCTIONS BEFORE USING (THIS APPLIANCE)**

**DANGER** - To reduce the risk of electrocution:

1. Always charge the product in a safe place away from water.
2. Do not reach for the product that has fallen into the water while charging. Unplug immediately.
3. Never charge with a damaged power cord or plug.
4. Make sure the charging port is dry enough before connecting it to the power supply.

**WARNING** - To reduce the risk of burns,

electrocution, fire, or injury to persons.

1. A product should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces.
6. Never use while sleeping or drowsy.
7. Never drop or insert any object into any opening or hose.
8. Do not expose the appliance or battery to excessive temperatures.
9. This appliance can be used by children aged from 12 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
10. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
11. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
12. The appliance must only be supplied at safety extra low

13. The power supply with output power not exceeding 15W in normal and abnormal working conditions must be used.
14. Be aware of the risk of terminals of the battery-operated appliance or battery being short circuited by metal objects.
15. If the battery is damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
16. This appliance is just for household use.

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