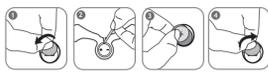
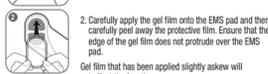
3.1 Inserting the battery

- Turn the battery compartment lid anti-clockwise using a coin.
- Use a pointed object to loosen the lid.
- Place the battery in the battery compartment so that the positive terminal (+) faces upwards.
- Close the lid by turning it clockwise.

3.2 Attaching the gel film

- Carefully remove the blue protective film.
2. Carefully apply the gel film onto the EMS pad and then carefully peel away the protective film. Ensure that the edge of the gel film does not get affected by the EMS pad.

Gel film that has been applied slightly askew will not affect the function.

⚠ WARNING!

- Make sure that you attach the gel film to the device on the side where the blue protective film was, otherwise the skin could become irritated.

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- Give the battery compartment lid anti-clockwise with the aid of a coin.
- Use a pointed object to loosen the lid.
- Place the battery in the battery compartment so that the positive terminal (+) faces upwards.
- Close the lid by turning it clockwise.

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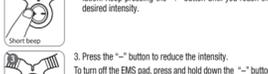
3.3 Usage

Stick the EMS pad onto the part of the body to be treated (see image).

- Press and hold down the "+" button until you hear a long beep.
- Press the "+" button again to start the level 1 stimulation. Keep pressing the "+" button until you reach the desired intensity.
- Press the "-" button to reduce the intensity.

- Turn off the EMS pad press and hold down the "-" button until you hear a long beep. The device switches off automatically after 20 minutes.
- Re-apply the protective film to the gel film after use.

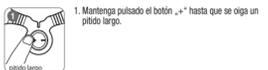
For optimal results we recommend a 4-week training plan. Please note the following points:

- 1-2 days rest.
- Increase the intensity level on a daily basis.
- A higher intensity level can lead to dizziness. However, the stimulation should never be uncomfortable.
- For optimum stimulation we recommend that the gel film is replaced when batteries are replaced.
- Do not use the EMS pad on the same area more than once per day.
- In the event of continuous burning, stop using the device. Clean the gel film and try to restart the program. If the continuous burning persists, replace the gel film. If the adhesive capability of the gel film decreases, replace immediately. Do not use the EMS pad again until the gel film has been replaced. Otherwise the unequal distribution of the gel film could lead to skin lesions.
- In the event of unpleasant current peaks, switch off the EMS pad, check the electrodes and then re-apply.
- The "Smart Contact Detection" automatically triggers an audible alarm if 2 electrodes come into contact with each other.

! Important

- The device will only work if there is sufficient contact with the skin. Stimulation pulses will not be transmitted if there is insufficient contact.
- Switch off the EMS pad before removing it.
- Replace the battery if it is low (insufficient stimulation, no beeping noises).
- For optimum stimulation we recommend that the gel film is replaced when batteries are replaced.
- Do not use the EMS pad on the same area more than once per day.
- In the event of continuous burning, stop using the device. Clean the gel film and try to restart the program. If the continuous beeping persists, replace the gel film. If the adhesive capability of the gel film decreases, replace immediately. Do not use the EMS pad again until the gel film has been replaced. Otherwise the unequal distribution of the gel film could lead to skin lesions.
- In the event of unpleasant current peaks, switch off the EMS pad, check the electrodes and then re-apply.
- The "Smart Contact Detection" automatically triggers an audible alarm if 2 electrodes come into contact with each other.

We reserve the right to make technical changes to improve and develop the product.

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On-time (s) = switch-on time in seconds (contraction).
Off-time (s) = switch-off time cycle in seconds (relaxation).

* There are 3 individual rounds in cycle 3.
** The off-time is 10 seconds for every 3rd round in cycle 3.

Note:The product components do not require any maintenance.

5. Disposal

The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronic retailers. You are legally required to dispose of the batteries. The codes below are printed on batteries containing harmful substances: Pb = Battery contains lead, Cd = Battery contains cadmium, Hg = Battery contains mercury.

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). See instructions for use. For more information, please contact the local authorities responsible for waste disposal.

6. Technical specifications

Typ		EM 20 EMS-Pad Neck
Output wave form	Biphasic rectangular pulse	
Pulse length	400 µs ± 10%	
Pulse frequency	4-6-50 Hz ± 10%	
Output voltage	max. 70 Vpp ± 10% (1000 Ohm)	
Output current	max. 70 mApp ± 10% (1000 Ohm)	
Voltage supply	3V CR2032 battery	
Treatment time	20 minutes, the device then switches off automatically	
Intensity	Adjustable from 0 to 15	
Operating con-	10°C – 40°C (50°F – 104°F) with a relative humidity of 30–85%	
Storage conditions	-10°C – 50°C (14°F – 122°F) with a relative humidity of 10–95%	
Dimensions	278 x 87 x 13 mm	
Weight	42 g (without battery), 45 g (with battery)	

- Note: If the device is not used according to the instructions specified, perfect functioning cannot be guaranteed!

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