Fairywill

   @FairywillCare


 www.fairywill.com

 support@fairywill.com

Thanks for choosing Fairywill Oral Irrigator. Before operating the product, please read the manual carefully and keep it for future reference.

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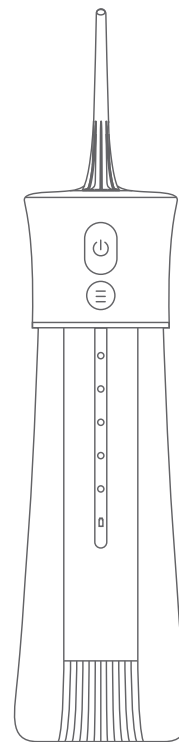Oral Irrigator

IMPORTANT SAFEGUARDS

All information in this section is provided to instruct you on the correct and safe operation of this device so that you can avoid injury to yourself and others as well as damage to property.


READ ALL INSTRUCTIONS BEFORE USING (THIS APPLIANCE)

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. This appliance contains batteries that are non-replaceable. When the battery is at the end of life, the appliance shall be properly disposed of.
4. Do not expose the appliance or battery to excessive temperatures.
5. Be aware of the risk of terminals of the battery-operated appliance or battery being short circuited by metal objects.
6. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
7. The power supply with output power not exceeding 15W in normal and abnormal working conditions must be used.
8. Toddlers, infants, people who are unable to operate the device should not use it by themselves. Please also keep it away and out of reach from children.
9. Never try to alter, disassemble or repair the device.
10. Do not fill the tank with water warmer than 40 °C.
11. If the external flexible cable or cord of this appliance is damaged, it shall be replaced by the manufacturer, a service agent or a similarly qualified person in order to avoid any hazard.
12. Keep the cord away from heated surfaces.
13. Please do not use it if you lose the nozzle.
14. Do not use for any purpose other than oral cleaning.
15. Keep nozzles away from young children. A nozzle is small enough to swallow and can be lodged in the throat.
16. Certain types of mouthwash may damage the equipment by cracking the housing or water reservoir. Therefore, do not use mouthwash in the water reservoir.
17. People who experience tooth or gum pain may not be able to use this device and should consult their dentist before use.
18. Do not press the nozzle too firmly against your teeth or gums, it will hurt your gums if you do so.

19. Please use an adapter so that the voltage is 5V $\overline{=}$ 1000mA during charge, otherwise the battery could be damaged by overcharging.

SAVE THESE INSTRUCTIONS

1. Designs:

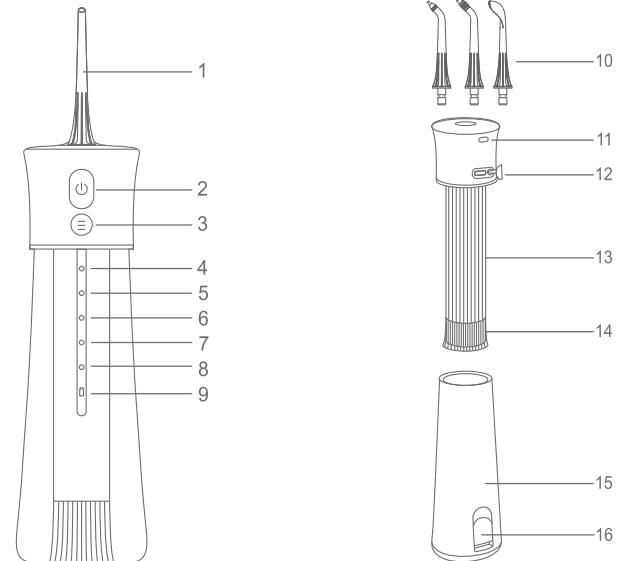- | | |
|-------------------|---|
| ① Standard Nozzle | ② On/Off Button |
| ③ Mode Button | ④ Soft |
| ⑤ Normal | ⑥ Potent |
| ⑦ Pulse | ⑧ Point |
| ⑨ Power Indicator | ⑩ Periodontal/Orthodontic/Tongue Cleaning Nozzles |
| ⑪ Unlock Button | ⑫ Charging Port (Charger Base) |
| ⑬ Flosser Panel | ⑭ Filter Part |
| ⑮ Water Tank | ⑯ Water Tank Lid |


Note:

There may be some liquid remaining inside the main device before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.


2. How To Use:

2.1 Charging the Device (As Fig. 1)

- Please cover the lid of the tank tightly before using the oral irrigator to avoid water leakage.
- Take the charger base cover down and insert the device socket (make sure no liquid is remaining inside the device plug).
- Insert the charging adapter securely into the power outlet.
- The “” indicator blinks when the device is charging. When it is fully charged, the white “” indicator will stay on for 2 minutes, and it will turn off automatically after that.
- Put the cover on the charger base again after charging is finished.


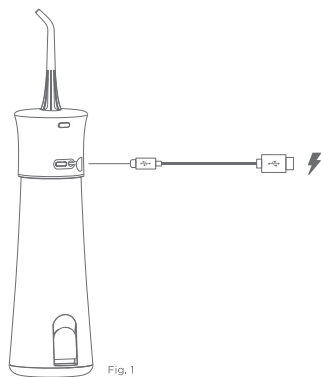
Note:


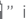
The “” indicator will flash in red when the power is lower than 10%; when the power is lower than 3%, the red “” indicator will flash for 5 times, and then it will turn off automatically, the device will also be powered off and is not able to work. It takes 4 hours to fully charge the device.

2.2 Use the Oral Irrigator

- Attach the nozzle onto the unit. (Fig. 2)
- Open the cap of the water tank to fill the tank with water which should not exceed 40°C, then close the cap.(Fig. 3)

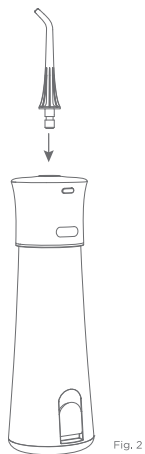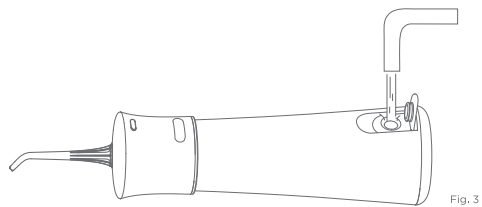- Or, slide the water tank down off the unit, fill the tank with water which should not exceed 40°C, then attach the tank on the unit. (Fig. 4)
- Press the “” button to change the working mode after the irrigator turns on, the default mode keeps the same as last shutoff in that situation. (Fig. 5)

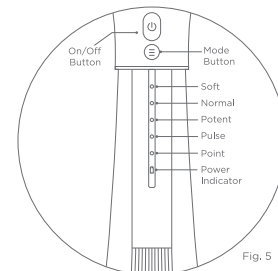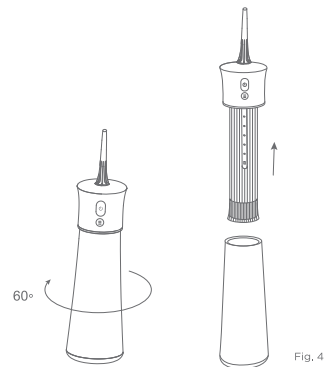Mode Introductions:

Soft: for gentle cleaning of sensitive gums, the water pressure will be 24-47psi.

Normal: for normal cleaning, the water pressure will be 40-71psi.


Potent: for people who have high requirements of teeth and gum cleaning, the water pressure will be 55-95psi.

Pulse: for massaging the gums, the water pressure will be around 24-47psi.


Point: for cleaning specific area of the teeth or gum, the water pressure will be 55-95psi, water will come out after pressing the button per time, but water will not come out when pressing the button for long time.

- Hold the unit in an upright position with the nozzle tip pointing toward your teeth.
- Keep your mouth opened slightly so that the water can flow out.
- Press the “” button shortly to turn on the device, it will start at the previous mode, you can also turn it off by pressing the “”.(Fig. 6)

- 1) Direct the water jet flow at a 90-degree angle to the teeth and gums.
- 2) Direct the water jet slowly along the teeth.
- 3) Direct the water jet along the gum line with the nozzle tip at a 90 degree angle to the gum tissue.
- 4) Cleaning between the teeth and around any braces, crowns or bridgeworks thoroughly.

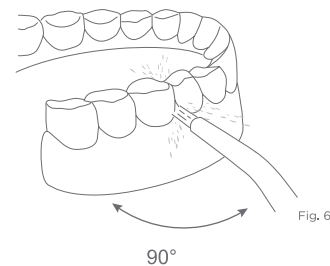- Firmly hold the base then, pull the nozzle out away from the base. (Fig. 7)
- Once irrigating is complete, press the “⏻” button again to turn off the device.

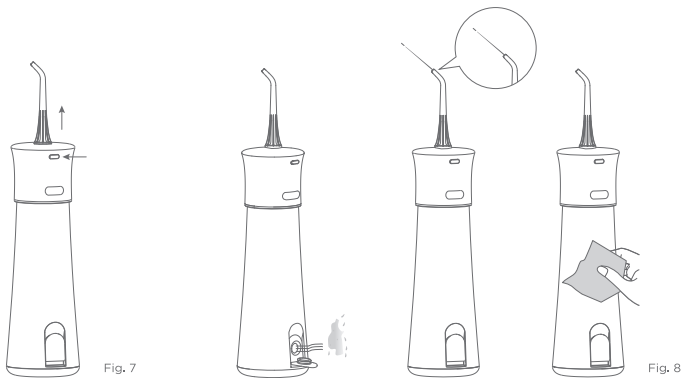2.3 After Using the Appliance (Fig.8)

- Open the cap of the water tank and dispose of the water in the tank.
- Press the “⏻” key to turn on the unit, and drain the water inside the device.
- Press the “⏻” key again to turn off the unit, and close the water tank cap.
- Wipe the device with a dry cloth.

Note:

To prevent bacteria from propagating, make sure there is no remaining water inside the unit. If you're not going to use this appliance for a long time, be sure to wipe and dry the appliance as soon as possible before storing it.

3. Specifications:

Working temperature: 0-50°C

Working voltage: DC3.7V

Input voltage: DC5.0V

Rated power: 6.0-7.0W

Waterproof rating: IPX7

4. Warranty:

- 1-year warranty is offered, starting from the first day of purchase. You may contact Fairywill.
- During the warranty period, for any defect that may arise from materials or product process, we will replace the parts or the whole product according to certain circumstances.
- When you apply for free repair, you must submit the product and proof that can prove the date of purchase to Fairywill.

- The following conditions are not included in the warranty:
 - Damages caused by improper use, for example, wrong use of power supply, plugged into an improper socket, artificial damage etc.
 - Ageing and normal wear and tear from use.
 - Problems occur due to its use for other purposes, or improper use or use components not supplied by Fairywill.
- This warranty does not affect your statutory rights.

5. Contact:

Fairywill

www.fairywill.com

support@fairywill.com

Made in China

Distributed by:

Zhuhai Zonghengyoupin Kuajingdianshang Youxiangongsi,
Zhuhai, Guangdong, 519031, CN

Recycling & Environmentally Friendly Disposal

There is a built-in rechargeable lithium battery in the product. When disposing of the appliance, do not mix it with general household waste. It should be sent to the battery recycling center for disposal in accordance with local disposal management regulations.

In order to protect the environment, after the expiration of the products service life, it should not be discarded in the same way as domestic garbage; please recycle at a professional recycling center. For further information on the recycling of this product, please consult your local municipality, processing service center, or the store where you purchased the product. You can help protect the environment! Please remember to respect the local regulations, and hand in the non-working electrical equipment to an appropriate waste disposal center.

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